

Product data sheet





Information


Designation Reference

White foldaway shower seat - Large size 047632

Product advantages

- + Maximum safety: with automatic stand.
- + Compact in raised position.
- + Easy fixing
- + Corrosion resistant warranty.

Description

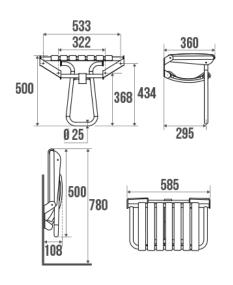

Large size. 7 white polypropylene slats. White epoxy-coated aluminium tube. Enables users to sit down safely and comfortably while showering.

Standard height 500 mm.

Technical characteristics

Material	Aluminium & Polypropylene
Finish	Epoxy-coated
Color	White
Dimensions	360 x 580 x 500 mm
Diameter Ø	25 mm
Gross weight	3.600 kg
Net weight	3.000 kg

PELLET ASC - ZI DE CHAPOTIN - 240 RUE DES FRÈRES VOISIN - CS 270048 - 69967 CHAPONNAY CEDEX - FRANCE - TEL : +33 (0)4 78 96 82 20 - FAX : +33 (0)4 78 96 70 52

© PELLET ASC – Non contractual document. Please read our General Sales Conditions. Dimensions in mm.

All dimensions and specifications are given for information purposes only and can be changed without notice in order to improve our products.



Packaging CE Marking : safe load Laboratory tests Guarantee Gencod Custom code Technical sheet Spare parts Cardboard box 150 kg 225 kg 10 years 3563390476321 7615200000 <u>A13600811_IND_15</u> Ref. 047636: pads for shower seat. Ref. 047638: White seat slat. Ref. 047639: Support slat.

Fixing and maintenance advices

Careful consideration and selection of the appropriate wall fixings according to the wall-mounted support is required to ensure the long life of your accessory. The use of stainless steel screws (Ø 4 mm) is highly recommended to prevent corrosion in damp environment. Before drilling, physically hold the product against the wall at desired position and mark the hole locations. All the fixing points must be respected. After installation and before use, make sure that all fixing screws are tightened securely. The shower seat provides support and added safety if properly fitted in accordance with our technical instructions. For use in seated position. To guarantee the user safety, DO NOT rise up on the seat. When the seat is raised, DO NOT use it as a grab bar. To avoid any risks of pinching out, DO NOT put your fingers in the mobile parts. The product does not contains any dangerous substance for the health. We disclaim all liability for any personal injury resulting from inappropriate installation or use of this seat.

Maintenance

IMPORTANT: To ensure the sustainability of the product and user safety, we recommend using a solution made of mild soap for cleaning. CAUTION: DO NOT USE cleaning products made of: etones, hydrocarbons, acids, bases, esters and ethers. After use of a cleaning product, make a rinsing in the clear water.

Confort of use

Function		Characteristics	Characteristics	
Rest	3/3	Horizontal	3/3	
		Vertical	3/3	
		Corrosion resistant	3/3	
Type of Use		Location of Use		
Shower	3/3	Medical facilities	3/3	
Cloakroom/Changeroom 3 /	3/3	ERP/IOP	3/3	
		Hotels	2/3	
		Private homes	3/3	
		Workplaces	3/3	

Fiche produit mise à jour le : 11/01/2019